SwiSto3

**NEU
NEW**

IONTOPHORESE HEIM- & PRAXISGERÄT

IONTOPHORESIS SET FOR THE TREATMENT OF HYPERHIDROSIS

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The elevated tendency to perspire on the feet, hands and underarms is often very uncomfortable for those affected. Iontophoresis has proven to be an effective, cost efficient in the long-term, safe and easy-to-use form of therapy for this situation. **Compared to injections or operative therapy, Iontophoresis therapy is pain-free. Medications can cause side effects.**

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SwiSto3 continuous and pulsed current therapy with one machine for the treatment of hand, foot and underarm perspiration.

KaWe's 3rd generation, the latest technology and manufacturing experience for more than 20 years.

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The new **SwiSto3**, with its modern design is extremely user-friendly and provides you with two possible treatment options for the treatment of hyperhidrosis:

- **therapy with continuous current or**
- **therapy with pulsed current,**

which is especially suitable for sensitive patients or children and youths. The new **SwiSto3** does not have language software. The user is simply guided through the menu and confirms the desired setting by pressing the SET button. The graphic LCD module illuminates shows all of the setting in a clear and organized manner. The user can change the current intensity unaided during the treatment, so that neither a second person nor a foot switch are necessary. The automatic current control function prevents the user of this equipment from feeling an electric shock (electric fence effect) when removing their hands and feet from the water trays during treatment. The **SwiSto3** is equipped with an automatic shutoff. If the device is not switched off after a treatment is finished, it switches itself off automatically after 2 minutes.

Scope of delivery

A charging plug is included and optionally, a UK-style plug can also be ordered. The **SwiSto3** can therefore be used in many different countries and due to its very light weight is an ideal travel companion. The Li-ion rechargeable battery is able to power up to 10, 20-minute treatments at full power on one charge.

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The plastic water trays and the electrodes are to be cleaned thoroughly after each use with disinfecting agent. For hygienic reasons, the mesh pads and sponge pockets are only to be used by one and the same patient. If possible, the use of skin crèmes and oils should be avoided during treatments.

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Features:

- no language software – can be used globally
- back-lit graphical LCD
- modern design, light-weight
- simple operation and menu navigation
- choice of continuous or pulsed current
- automatic current control, no feeling of electric shock when hands and feet are removed from the water trays
- treatment of the hands without foot switch possible
- membrane keypad, fully-automatic operation
- timer function
- automatic shut-off, sleep mode
- Li-ion rechargeable battery with charging function, charging plug
- per charge: up to 10, 20-minute treatments
- easy to disinfect
- optional UK-style charging plug available
- complies with all MPG requirements and standards
- medication-free therapy
- **4 year warranty** except for the Li-ion rechargeable battery

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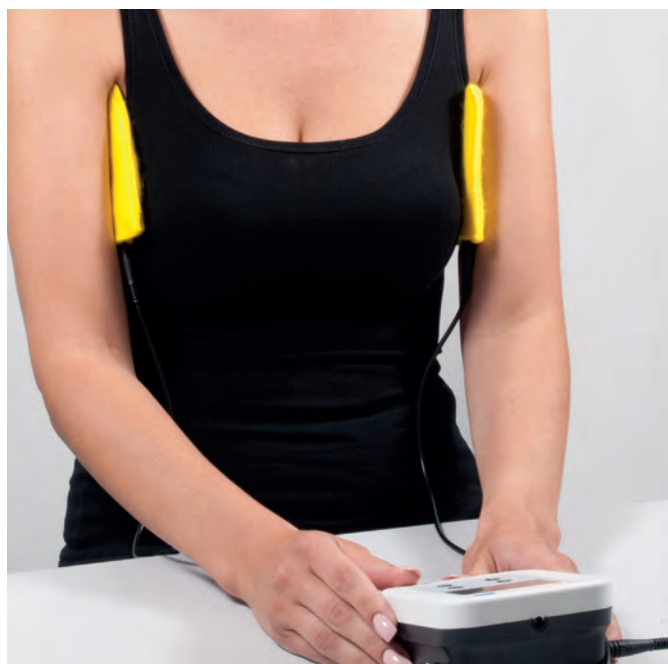
Use

Iontophoresis treatment for hyperhidrosis manum et pedum is a physical application. It was first described by Boumann and Grunewald-Lentzer in 1952. **Continuous or pulsed current** cause ion transport through the skin.

The plastic treatment trays are filled with approx. 3-4 cm of luke-warm water, such that the surfaces of the hands and feet are covered. Before filling the trays, the electrodes must be covered with the mesh pads in order to prevent direct contact with the skin. The sponge pockets used for underarm treatment are to be soaked thoroughly with water before use. The treatment should be performed at least 3 times per week with an average treatment time of 15 to 20 minutes until the desired healing results have been achieved. Usually, this takes place after approx. 10-15 treatments. Thereafter, maintenance treatments should be performed regularly (at least once a week) in order to maintain the successful therapy results. The current that flows through the electrode plates can, depending on application (hands, feet, underarms), be adjusted up to a maximum of 25 mA. We recommend, that the polarity be switched if the skin is thickly calloused as the anode is said to have a stronger empirical efficacy.

Additional healing results

In addition, scientific test have shown that this therapy can also lead to excellent treatment results for gram-negative foot infection, dyshidrotic eczema, pustulosis palmaris et plantaris, acrocyanosis and osteoarthritis.

Please consult your physician if you have any of the following:

- Cardiac arrhythmia
- Pregnancy
- Pacemakers or metal implants
- Metal-containing intrauterine devices (contraceptive coil)
- Major skin defects that cannot be covered with Vaseline
- Insensitivity to painful stimuli

SwiSto3

GLEICH- UND PULSSTROM | CONTINUOUS AND PULSED CURRENT

2in1

SwiSto3 – Iontophorese-Set SwiSto3 – continuous and pulsed current

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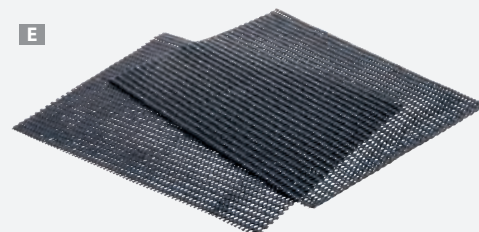
Iontophoresis device as replacement device
with plate electrodes for the treatment
of underarm perspiration

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
made of stainless-steel,
incl. sponge pockets,
incl. cable

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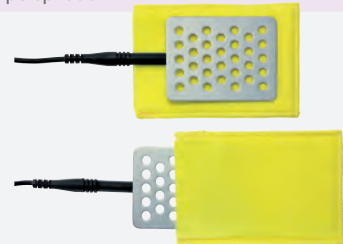for plate electrodes,
For the treatment of underarm perspiration

Maße | Dimensions
ca. | approx. 11 x 9 cm

(VE = 2 St.) | (PU = 2 items)

REF 05.19090.001 (68580)

133

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Behandlungswanne Water trays

Maße | Dimensions
ca. | approx. 37 x 31 x 7 cm

REF 05.19020.021

Elektroden-Netzauflagen Electrode mesh pads

Maße | Dimensions
ca. | approx. 25 x 20 cm

(VE = 2 St.) | (PU = 2 items)

REF 05.19050.001 (68531)

Elektroden-Platten Electrode plates

inkl. Kabel | incl. cable

Maße | Dimensions
ca. | approx. 22 x 15 cm

(VE = 2 St.) | (PU = 2 items)

REF 05.19040.021

Ladestecker | Charging plug

ca. | approx. 0,10 kg

240 V
EU
VERSION

240 V ~ 50-60 Hz, 400 mA

REF 05.19170.002