



Rapid Rhino™

Epistaxis Solutions



Designed to Create:

- Fast, single treatment solution
- Low profile and self-lubricating
- Inflatable Tamponade

Rapid Rhino

Offering a Full Range of Epistaxis Solutions


Ultra low profile and self-lubricating hydrocolloid fabric provides for unsurpassed¹ ease of insertion and removal.




CMC - infused exterior fabric wets into a self-lubricating gel for optimal patient comfort from insertion to removal

Key Features:

- Easy to insert and remove
- Quickly conforms to nasal anatomy
- Provides gentle and even compression to areas of epistaxis
- Designed to prevent adherence to tissue or blood clots upon removal
- Available in a wide range of sizes

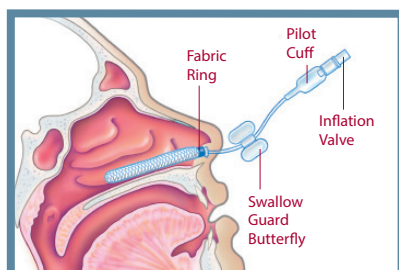

Rapid Rhino Nasal Pac
(Uninflated 5.5cm Anterior
Shown)

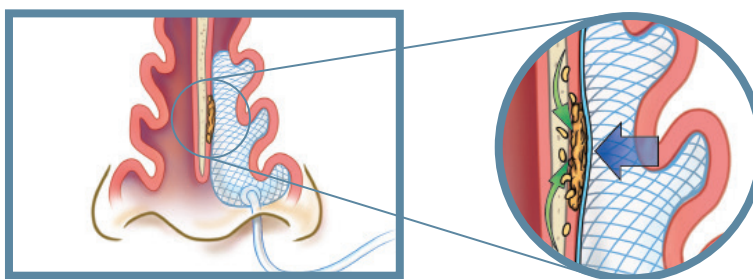


Rapid Rhino with
Carboxymethylcellulose (CMC)


A Comparison

A comparison of Rapid Rhino nasal tamponades to Rhino Rocket™² for the treatment of epistaxis in the emergency department¹: Rapid Rhino epistaxis products provide greater patient comfort at the insertion and at removal.


Ease of Insertion

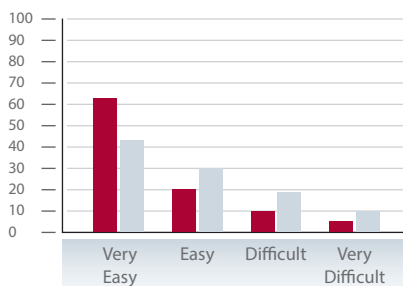

Ease of Removal

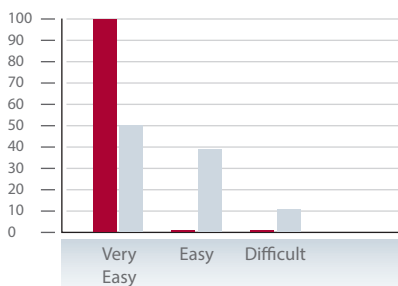

Mean Pain Scores

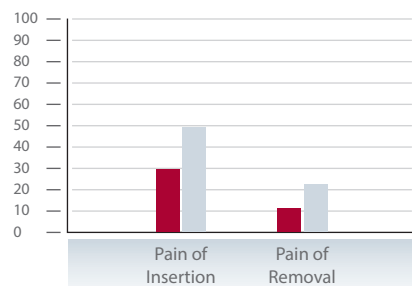


1. Singer AJ, Blanda M, Cronin K, LoGiudice-Khwaja M, Gulla J, Bradshaw J, Katz. A Comparison of nasal tampons for the treatment of epistaxis in the emergency department: a randomized controlled trial. *Ann Emerg Med.* 2005; 45:134–139.

■ Rapid Rhino
■ Rhino Rocket


Rapid Rhino Product Usage Directions




- Soak in sterile water for a FULL 30 seconds.

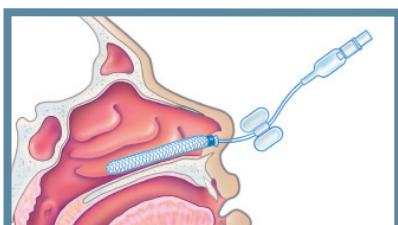


- Insert along superior aspect of the hard palate until the blue indicator is past the nares.

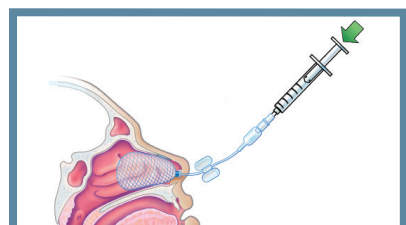


- Using a 20ml syringe, inflate the Rapid Rhino device with AIR only. Monitor the pilot cuff for direct tactile feedback.
- Stop inflation when the pilot cuff becomes rounded and feels firm when squeezed.

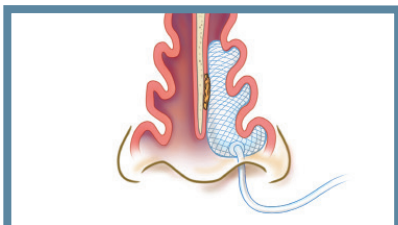

- Inflate the cuff to provide a gentle, low-pressure tamponade delivering the CMC fabric directly to the bleed site.

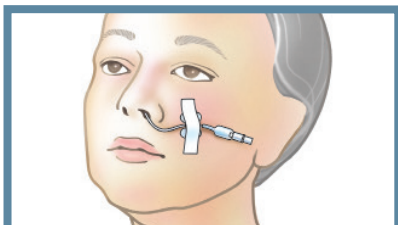

- Reassess after 15-20 minutes.
- Reinflate to ensure proper pressure (if necessary) and tape to patient's cheek away from the upper lip.

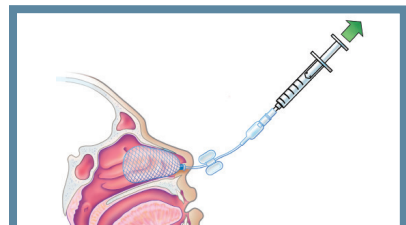

- Removal should occur 24-72 hours after treatment.

Ordering Information

Part Number	Description
RR450	4.5cm Inflatable Unilateral without airway
RR530	5.5cm Non-Inflatable Rapid-Pac™
RR550	5.5cm Inflatable Unilateral without airway
RR551	5.5cm Inflatable Unilateral with airway
RR750	7.5cm Inflatable Unilateral without airway
RR751	7.5cm Inflatable Unilateral with airway
RR900	9cm Inflatable Dual Balloon

To order and for more information, please contact your ArthroCare representative



Rapid Rhino 450

- 4.5cm anterior inflatable
- For use in patients with small nasal anatomy



Rapid Rhino 550

- 5.5cm inflatable epistaxis device



Rapid Rhino 551

- 5.5cm anterior nasal tamponade that contains an airway



Rapid Rhino 750

- 7.5cm inflatable epistaxis device



Rapid Rhino 900

- 9cm tamponade device for posterior epistaxis
- Contains 2 balloons to control intranasal hemorrhage originating from the sphenopalatine artery



Rapid Rhino 751

- 7.5cm anterior/posterior device that contains an airway



Rapid Rhino 530

- Rapid-Pac™ is 5.5cm
- Offers a non-inflatable option for minor anterior nosebleeds
- Treats less severe bleeds in the ER or urgent care



arthrocare.com/international

Corporate Headquarters

7000 West William Cannon Drive
Austin, TX 78735
Main Phone: +1.512.391.3900
Order Entry: +1.800.797.6520
Order Entry Fax: +1.888.994.2782

EMEA Headquarters

St. James Business Park
9, Grimbald Crag Court
Knaresborough
North Yorkshire HG5 8QB
Phone: +44 (0) 1423 888806
Fax: +44 (0) 1423 888816

Asia Pacific Headquarters

17th Floor, Tower B
Ping An International Financial Centre
No. 1-3, Xin Yuan Nan Road
Chao Yang District Beijing 100027
China
Phone: +852 3793 3430 32

Authorised Representative

Skeppsbron 2
111 30 Stockholm
Sweden
Phone: +46 8 546 172 00
Fax: +46 8 546 172 39

²Rhino Rocket is a registered trademark of Shippert Medical Technologies Corp. and is used for reference purposes only.

Not For Distribution in the United States.

©2013 ArthroCare Corporation. All rights reserved. All trademarks, service marks and logos listed herein are registered and unregistered trademarks of ArthroCare Corporation, its affiliates or subsidiaries or a third party who has licensed its trademarks to ArthroCare.