



PRODUCTS FOR SALIVARY DUCT ACCESS

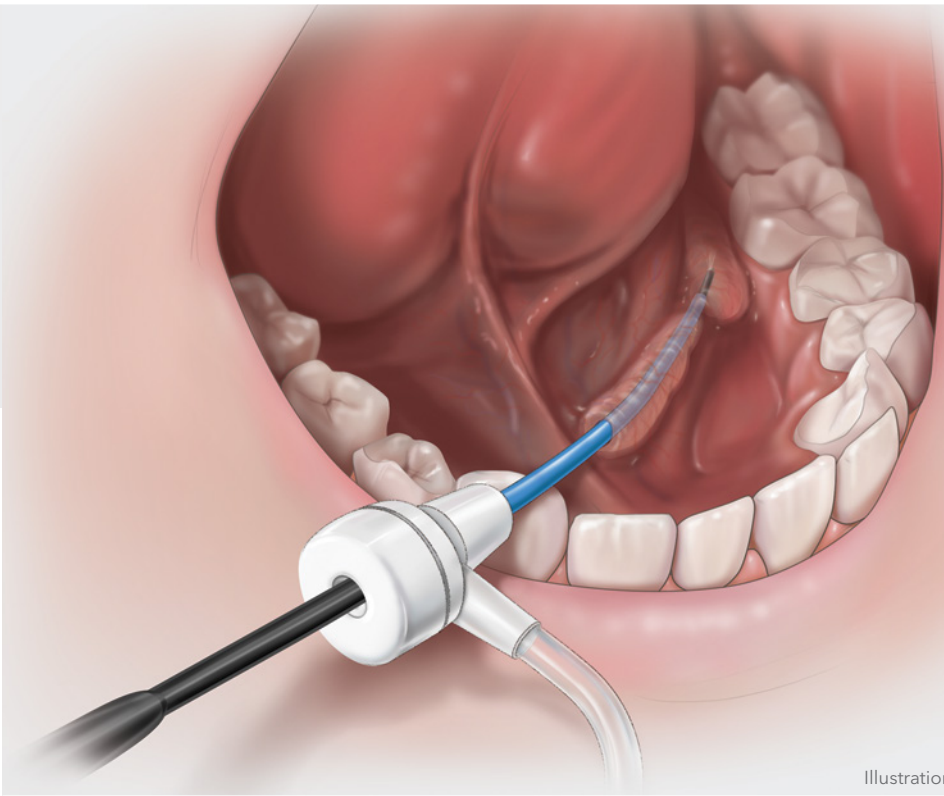


Illustration by Lisa Clark

Salivary Duct Access

Gain and maintain access to the salivary duct while reducing trauma and maximizing efficiency.

Sialendoscopy represents an exciting **minimally invasive** alternative to a standard open surgical procedure to treat salivary duct disorders. Sialendoscopy offers a variety of advantages to the patient, including minimal recovery time, no incision or subsequent scarring, minimized potential for nerve damage, preservation of the salivary gland¹ and early return to a normal diet. Sialendoscopy allows **stones to be removed in the office** with a low risk of major complications to the patient.

Cook Medical continually collaborates with physicians to improve procedures. Cook has partnered with leading physicians to create salivary duct access products that are designed to maintain access throughout the procedure and that work together for greater control. This **control** can help reduce trauma for the patient and increase efficiency with every step.

“Cook’s salivary duct access tools allow me to gain and maintain access to the duct, which gives me much greater control throughout the entire procedure.”

Jack Kolenda, MD
Chief, Department of Otolaryngology
Oakville-Trafalgar Memorial Hospital
Oakville, Ontario


GAIN AND MAINTAIN ACCESS

1

WIRE GUIDE

Gain and maintain access with a wire guide.

Sialendoscopy starts with accessing the salivary duct through the papilla. A flexible Cook soft-tip wire guide can be easily maneuvered into the duct. This guide allows you to easily gain and maintain ductal access throughout the procedure, which reduces trauma and maximizes efficiency.

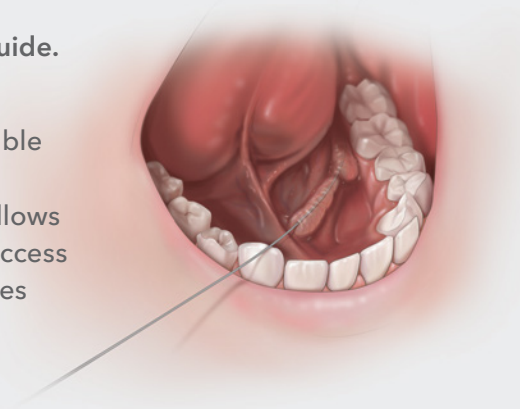


2

DILATOR SET

Minimize trauma with sequential dilation.

Once you establish access using the Cook soft-tip wire guide, you can introduce a series of flexible dilators over the wire in order to easily expand the papilla and to prepare the salivary duct for the introduction of procedural instruments.

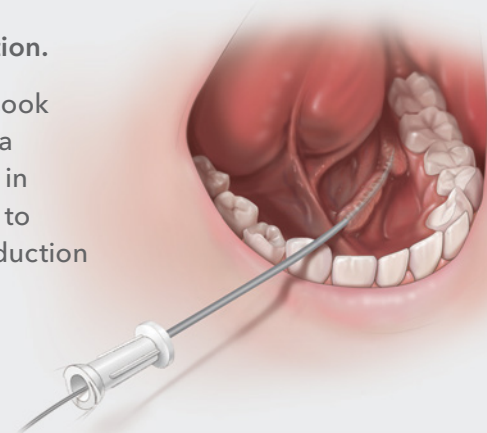


3

KOLENDA INTRODUCER SHEATH

Create a true working channel.

After you expand the papilla using the dilator set, you can create an open working channel into the salivary duct by passing the sheath over the wire. The sheath also serves to protect the ductal wall, which allows you to easily insert and remove procedural instruments.

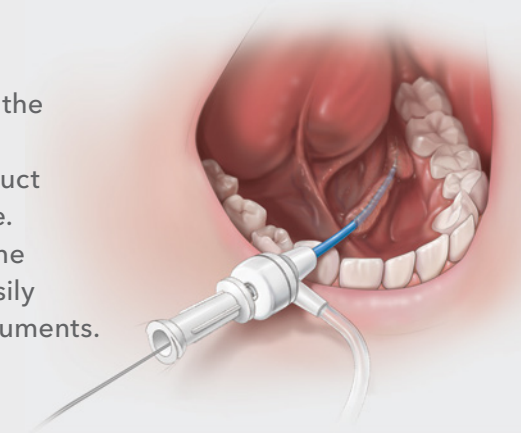

Illustrations by Lisa Clark

WIRE GUIDE

Order Number	Reference Part Number	Diameter inch	Length cm	Tip Configuration
G02255	C-SF-18-40	.018	40	straight
G02274	C-SF-15-50	.015	50	straight

ADAPTOR

G05301	PMLLA-VTA-L
--------	-------------

SALIVARY ACCESS DILATOR SET

Order Number	Reference Part Number	Dilator OD Fr/mm/inch
G24131	SDDS-18-20-HC	4.0/1.35/.053
		5.0/1.67/.066
		6.0/2.0/.079
		7.0/2.33/.092

KOLENDA SALIVARY ACCESS INTRODUCER SET

Order Number	Reference Part Number	Introducer ID mm/inch	Wire Guide Diameter inch	Sheath Length cm
G24132	SDAS-5.0-18	1.9/.074	.018	5.5
G24129	SDAS-6.0-18	2.2/.087	.018	5.5

Reference

1. Wallace E, Tuzin M, Hagan JP, et al. Management of giant sialoliths: review of the literature and preliminary experience with interventional sialendoscopy. *Laryngoscope*. 2010;120(10):1974-1978.

Resources

1. Chossegros C, Guyot L, Richard O, et al. A technical improvement in sialendoscopy to enter the salivary ducts. *Laryngoscope*. 2006;116(5):842-844.
2. Walvekar RR, Razfar A, Carrau RL, et al. Sialendoscopy and associated complications: a preliminary experience. *Laryngoscope*. 2008;118(5):776-779.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do

Distributors: +353 61239240, ssc.distributors@cookmedical.com

Austria: +43 179567121, oe.orders@cookmedical.com

Belgium: +32 27001633, be.orders@cookmedical.com

Denmark: +45 38487607, da.orders@cookmedical.com

Finland: +358 972519996, fi.orders@cookmedical.com

France: +33 171230269, fr.orders@cookmedical.com

Germany: +49 6950072804, de.orders@cookmedical.com

Hungary: +36 17779199, hu.orders@cookmedical.com

Ireland: +353 61239252, ie.orders@cookmedical.com

Italy: +39 0269682853, it.orders@cookmedical.com

Netherlands: +31 202013367, nl.orders@cookmedical.com

Norway: +47 23162968, no.orders@cookmedical.com

Spain: +34 912702691, es.orders@cookmedical.com

Sweden: +46 858769468, se.orders@cookmedical.com

Switzerland - French: +41 448009609, fr.orders@cookmedical.com

Switzerland - Italian: +41 448009609, it.orders@cookmedical.com

Switzerland - German: +41 448009609, de.orders@cookmedical.com

United Kingdom: +44 2073654183, uk.orders@cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do

Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335

E-mail: orders@cookmedical.com

Australia:

Phone: +61 734346000, 1800777222, Fax: +61 734346001, 1800077283

E-mail: cau.custserv@cookmedical.com

www.cookmedical.com

